

USER MANUAL
EVOLI SCALE
EVSC-5HB

Read carefully before using this product

Important safety precautions when using electrical appliances, certain basic safety precautions must be followed, including the following:

1. This product must not be used by children or by other people without assistance or supervision if their physical, sensory or mental condition limits safe use of the appliance.
2. In order to ensure reliable readings and guarantee connectivity, please stand on the scales with bare feet. Please make sure that both the scales and your feet are completely dry before use.
3. Never place the scales on carpet or on uneven floor.
4. The scales are not waterproof, never immerse the scales in water.
5. To clean the surface, use a damp cloth or a glass cleaner. Do not use soap or other chemicals. Please keep the scales away from water, heat and extreme cold.
6. The scales are a high precision measuring device never jump or stop on the scale or disassemble. Please handle the scales carefully to avoid breakage.
7. The scales are for domestic use only. You should consult a doctor when undertaking any diet or exercise program.
8. Remove the batteries if the scales are not being used for a long period of time.
9. Please check the batteries if the scales malfunction. Change to new batteries if required.
10. Please remove insulating strip on the batteries before use.
11. If batteries leak, change them immediately as they may damage the scales

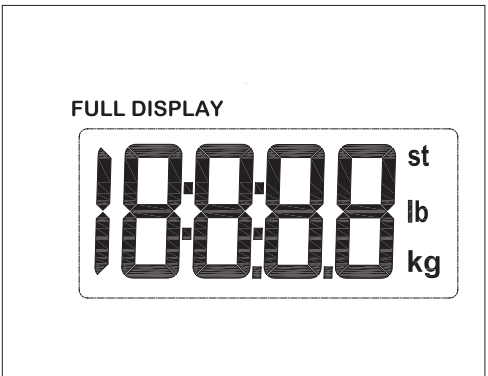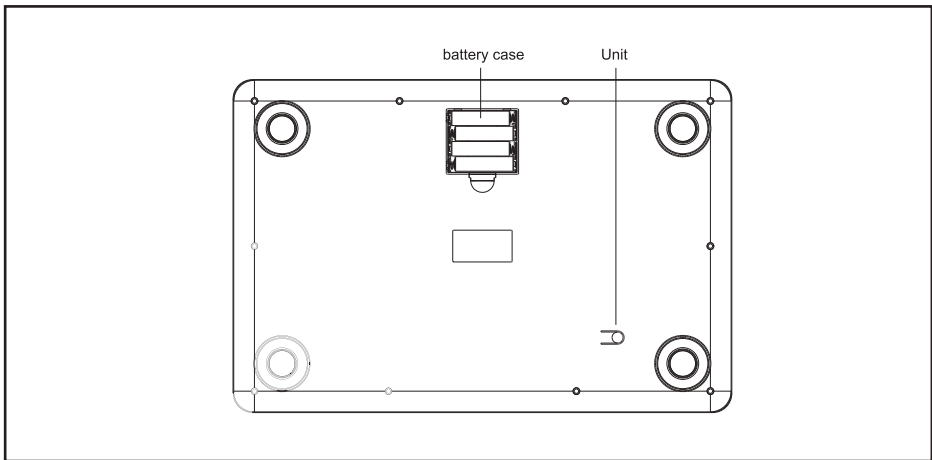FEATURES AND SPECIFICATIONS

1. auto on/off
2. high precision strain gauge sensors system
3. weight units: kg, lb, st:lb
4. maximum capacity: 150kg/330lb/24st
5. tempered glass platform – 361mm x 241 x 6mm
6. large LCD display with backlit easy to read – 110x 45mm
7. graduation: 0.1kg/0.2lb
8. overload indication: “Err”
9. low battery indication “Lo”
10. power supply: 4xAAA batteries
11. batteries included
12. one-year warranty

OPERATING INSTRUCTIONS

1. Remove all packaging before use.
2. Locate the 4 x AAA batteries and the battery holder at the back of the scales, open the cover, install the batteries into the holder and close the cover. Please ensure that the batteries are placed the correct way as indicated in the holder.
3. Press UNIT key to switch on and select unit measurement.
4. Place the scales on solid, flat ground, let the scales turn off automatically
5. Step on the scales gently, position your feet evenly near the center and keep your body still and erect until the scales stabilize your weight. The scales will flicker and lock on the LCD screen after your weight is steady, then automatically turn off after 10 seconds.
6. The scales will automatically turn off if no further use is required.

WORKING CONDITIONS

1. Operating temperature range: 0 - 40 °C.
2. Operating humidity range: ≤ 90% RH.
3. Storage temperature: -20 °C - 60 °C.

Meaning of crossed – out wheeled dustbin:

Do not dispose of electrical appliances as unsorted municipal waste, use separate collection facilities.

Contact you local government for information regarding the collection systems available.

If electrical appliances are disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being.

When replacing old appliances with new ones, the retailer is legally obligated to take back your old appliance for disposals at least free of charge.

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